

HG Series Totaldraft Stations



HG100 Totaldraft, Right Sink Approach



Totaldraft chamber

The HG Series Grossing Workstations are ventilated with Mopec's Total Downdraft System. This unique ventilating system moves air consistently and continuously across the entire work surface.

Stainless steel counter is fabricated with formed edges that are 1-1/4" thick with a marine edge around the perimeter. The HG100 and HG200 Inner chamber is equipped with a 14" x 12" x 9" deep stainless steel drain sink made of 14 gauge stainless steel. The work surface features self supporting perforated grid plates and is fabricated of 16 gauge stainless steel and perforated with 1/2" diameter holes on 1" centers. Grid plates are designed for easy cleaning with single flanges sealing off any unseen areas that could accumulate bacteria. Finger holds are included for easy handling. TIG welding on all seams and mating sections and ground smooth and regrained to match the #4 satin factory finish. All internal plumbing and electrical lines are factory installed and tested. Utilities hookup is minimal and can be accomplished in minutes.

Standard Features

Support Pedestal: All stainless steel construction. Fabricated of 18 gauge stainless steel, the support pedestal is complete with access panel. The cabinet is

fabricated for heavy-duty use. Material finish is #4 satin finish.

Hot/Cold Fixture: Hot/cold water fixture with swing spout and vacuum breaker.

Hand Spray: Soft spray head for easy positioning .

Electrical Receptacle: Electrical receptacle is provided. The receptacle is GFCI type with watertight cover.

Optional Features

Waste Disposal: BL800: Heavy-duty 1 HP commercial disposal, complete with on/off switch, water solenoid, and vacuum breaker plumbed and wired ready for use.

Hand Spray Rinse: Soft spray head for easy positioning.

Unit can be customized to fit your lab requirements with other additional options. Contact Customer Service with your specific needs.


HG100 and HG200 Totaldraft Stations



HG300 Totaldraft, Right Sink Approach

HG Series Totaldraft Stations

L = Length	Sink	Product No.
60" (152.4 cm)	Right	HG100
60" (152.4 cm)	Left	HG200
72" (182.9 cm)	Right	HG300
72" (182.9 cm)	Left	HG400

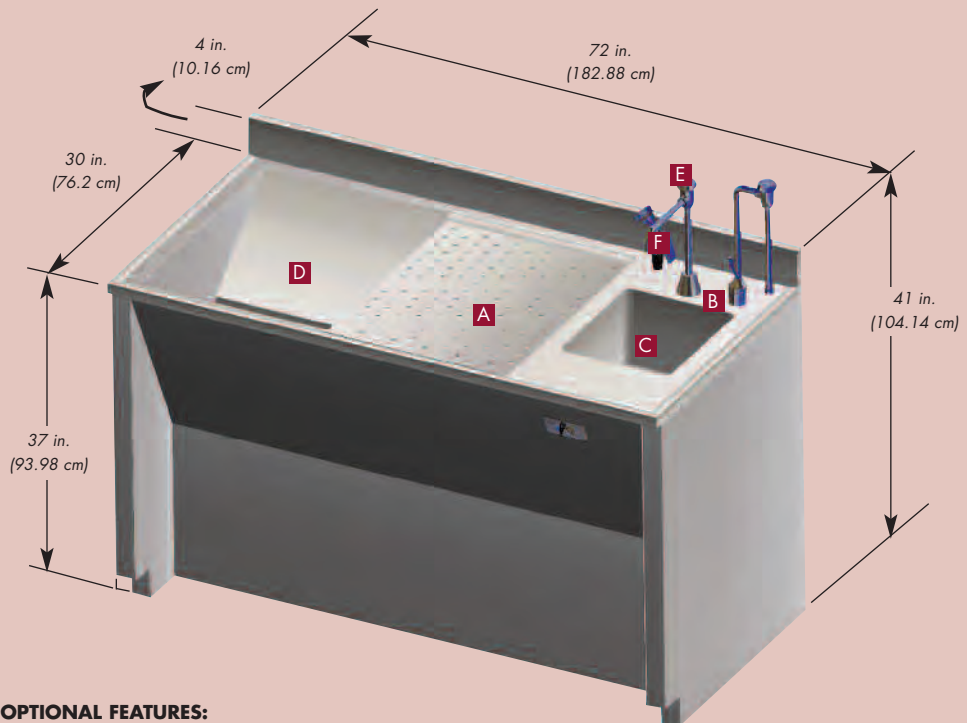

STANDARD FEATURES:

- A** Perforated Grid Plates
- B** Hot/Cold Mixing Valve
- C** Sink Well
- D** Totaldraft Ventilation Chamber
- E** Hot/cold fixture

OPTIONAL FEATURES:

- BL800 – 1 HP Disposal (Not Shown)
- F** Hand Spray Rinse

HG100 Shown with (1) one of (2) two Perforated Grid Plates.