

Fluid Collection



Vacuum Hydro-Aspirator

- High vacuum maintains a 26" Suction at 60 lbs. of water pressure
- Can easily disassemble for cleaning
- Built-In vacuum breaker eliminates the possibility of back siphonage
- Reverse flow allows for simple clearing of clogged suction line and provides fresh water spray 3-5 gallons per minute, 40-60 PSI
- Standard 3/4" (1.9 cm) hose thread for easy attachment to most sinks
- Brass construction with high luster chrome finish


Product No. BE013



Non-Clogging Postmortem Aspirator Extension Accessory

- Postmortem aspirator is used in conjunction with the Hydro Aspirator to evacuate fluids from the body
- Flat head allows complete drainage of cavities
- Can be used for Irrigation
- Brass construction with high luster chrome finish
- 10-1/2" (26.7 cm) Long

Product No. BE014

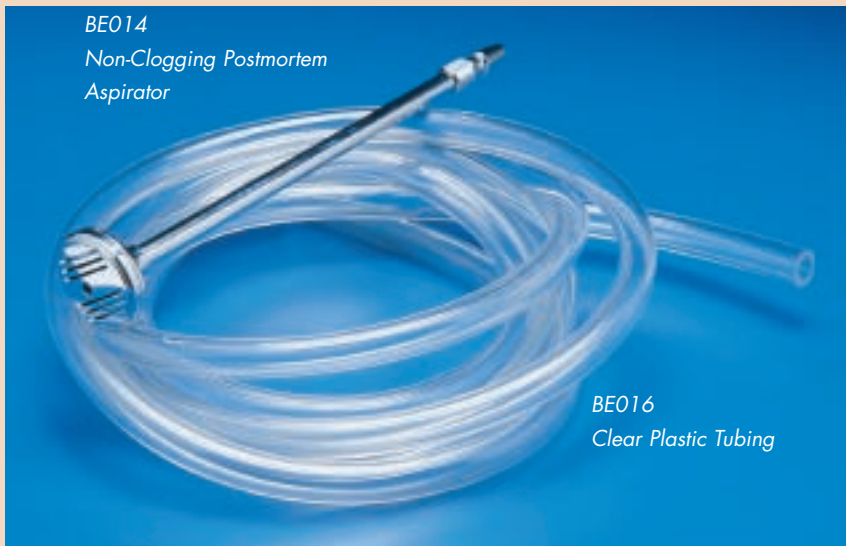

Clear Plastic Tubing - Thick Walled

- Works with most morgue plumbing fixtures
- Flexible
- 5/8" (1.5875 cm) Outside diameter
- 3/8" (.9525 cm) Inside diameter
- 1/8" (.3175 cm) Thick walls
- Sold in 10 ft. increments (25.4 cm) or greater

Product No. BE016



Hydro-Aspirator Extension

- The extension serves a dual purpose.
- Adds an additional vacuum breaker for added protection
- Simple attachment to the top of the Hydro Aspirator
- Brass construction with a high Luster chrome finish
- 5" (12.7 cm) long

Product No. BE015

The Mopec fluid collection can easily be added to any existing hydro-aspirator