


Laboratory Freezers/Refrigerators

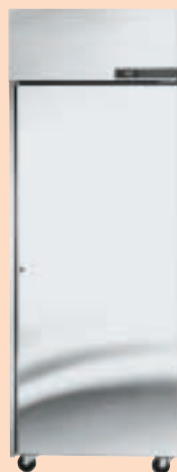

Laboratory Storage Unit with Digital Thermometer LED

All laboratory storage units come standard with digital thermometer LED (°C), complete perimeter anti-condensate door opening heater wires to prevent "freeze-out" and a door activated interior fluorescent light. Door mounted security lock. Includes adjustable temperature control range. Although all models come standard with 3 epoxy coated wire shelves, up to 5 additional shelves may be purchased.

Units feature a top-mounted, non-CFC refrigeration system (R404A) with an air-cooled condensing unit and automatic defrost. Adjustable temperature control range is also standard. Condensation is automatically evaporated in the collection pan.

Cabinet interior and exterior door and sides are fabricated of premium grade stainless steel. Galvanized exterior top and bottom. Not For Export.

1-Door Specifications

Dimensions:	81.5" H x 27.5" W x 34.6" D (207 cm x 69.9 cm x 87.6 cm)
Weight:	375 lbs. (170.1 kg)
Gross Cubage:	22.5 ft.3 (.64 m3)
Electrical:	115/60/1
Maximum Fuse Size:	15
NEMA Plug Config.:	5-15P
Total Amp Draw:	12
Refrigerator	Product No. KG100
Freezer	Product No. KG101

2-Door Specifications.

Dimensions:	81.5" H x 55" W x 34.6" D (207 cm x 139.7 cm x 87.6 cm)
Weight:	660 lbs. (299.3 kg)
Gross Cubage:	48.8 ft.3 (1.38 m3)
Electrical:	115/60/1
Maximum Fuse Size:	20
NEMA Plug Config.:	5-20P
Total Amp Draw:	14
Refrigerator	Product No. KG200
Freezer	Product No. KG201

3-Door Specifications

Dimensions:	82.75 H x 82.5" W x 34.6" D (210.2 cm x 209.6 cm x 87.6 cm)
Weight:	805 lbs. (365.1 kg)
Gross Cubage:	75 ft.3 (2.12 m3)
Electrical:	115/208-230/60/1
Maximum Fuse Size:	15
NEMA Plug Config.:	N/A*
Total Amp Draw:	11
*This model is available with a cord connection and casters. Consult Mopec.	
Refrigerator	Product No. KG300
Freezer	Product No. KG301