

# HG500 & HG600 Totaldraft Station

HG500 - Sitting 30" / HG600 - Standing 37"



*Our economical HG Series Grossing Stations can be customized to accommodate special requirements*

## HG500 & HG600 Downdraft Grossing Workstations

This unique downdraft table includes a multi-functional surface for specimen dissection. The large open work area has ample countertop to layout all dissection materials. This economical grossing station is all stainless steel construction and can be customized to accommodate special requirements.

### Standard Features

**Totaldraft System Design:** Totaldraft system will maintain a consistent airflow distribution over the full

working area. Duct stubs are included for connection to facility ventilation system.

**Gridplate:** Grid plate is designed for flush fit with table edge.

**Hot/Cold Water Fixture:** Swing spout faucet with built in vacuum breaker, has easily accessible wrist blade handles for hot and cold water control.

**Product No. HG500 - Sitting**

**Product No. HG600 - Standing**

## HG500 & HG600 Totaldraft Station

HG500 - Sitting 30" / HG600 - Standing 37"




### Optional Features

**Disposal:** BL800 1 HP commercial grade disposal

**Hand Spray Rinse:** Soft spray head for easy positioning.

Unit can be customized to fit your lab requirements with other additional options. Contact Customer Service with your specific needs.

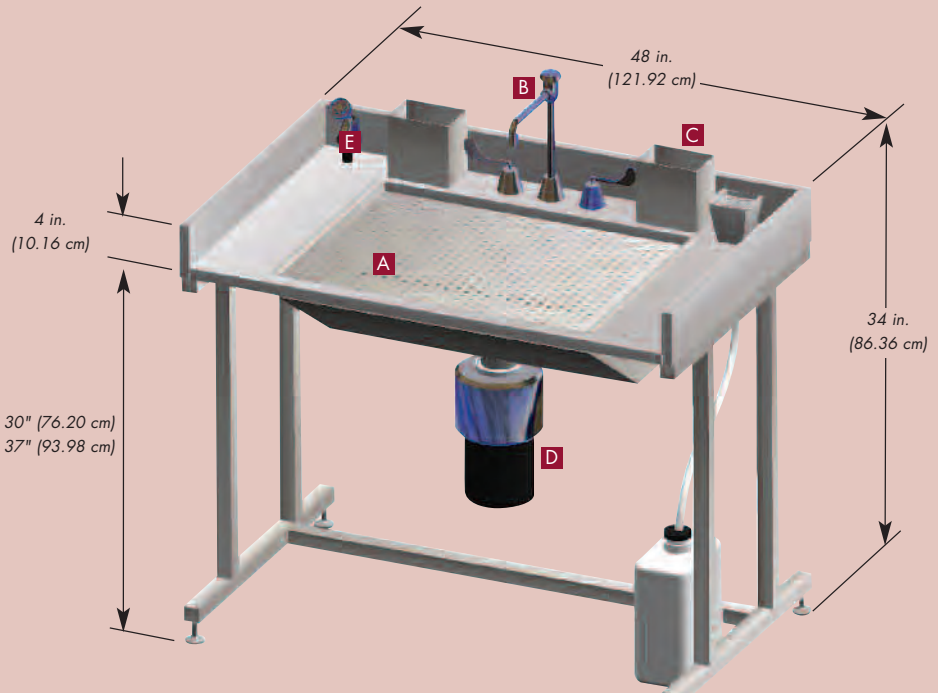

#### STANDARD FEATURES:

- A** Perforated Grid Plate
- B** Hot/Cold Water Fixture
- C** (2) Vent Stubs

HG500: 30" (76.20 cm)  
HG600: 37" (93.98 cm)

#### OPTIONAL FEATURES:

- D** BL800 1 HP Disposal
- E** Hand Spray Rinse