

# MB Series Optional Features

## Customized to Meet Any Laboratory Specifications

### MB001 Self Contained Exhaust System for MB100 and MB200

Dual blower air filtration system, complete with pilot lighted on/off switch, variable speed control on one blower only and filter caddy. Includes two potassium permanganate filters which allow for independent means of filtering the ventilated air, eliminating the need to tie into the building ventilation system. Requires (2) filters.

### MB003 Formalin Dispensing Carboy

2.5 gallon capacity Nalgene polyethylene dispensing container. Formalin not included.

### MB004 Dictation Equipment Stand

Angled surface for easy viewing of most dictation systems. Dictation stand is portable and can be placed anywhere. Stainless steel construction. Mouse not included.

### MB005 Self Contained Exhaust System for MB400 and MB600

Dual blower air filtration system, complete with pilot lighted on/off switch, variable speed control on one blower only and filter caddy. Includes three potassium permanganate filters. This system allows for an independent means of filtering the ventilated air, eliminating the need to tie into the building ventilation system.

### MB006 Leg Frame

Leg frame weldment designed to support the MB100 and MB200 grossing workstation allowing the grossing workstation to be free standing. Fabricated of 1-1/2" diameter heavy wall stainless steel tubing with adjustable leveling feet.

### MB007 Foot Pedal

Includes a dual foot pedal H/C water control valve and swing spout faucet for hands free control. Wrist blade handle controls are eliminated with this option.

### MB008 Flammable Storage Cabinet for MB400 only

Lower base cabinet designed to hold flammable liquids. Adjustable shelving included for easy and versatile storage, 22 gallon capacity. MB008 option is **not** compatible with options MB021 and/or MB023.

### MB009 Perimeter Rinse

Full perimeter rinse on three sides virtually eliminating the accumulation of debris and fluids in the dissection area. Control valve included for water flow adjustment.

### MB011 Safety Splash Shield

Lexan shield on flexible arm for easy positioning provides an additional shield between you and the specimen. The 11" x 12" contoured viewing area allows for a safe and large working area.

### MB012 Hand Spray Rinse

Soft spray head for easy and controlled rinsing of work surface. Chrome plated brass hand sprayer with 60" high-pressure hose for heavy-duty use. Control valve included for easy control of water pressure along with atmospheric vacuum breaker protection.

### MB013 Halogen Light Magnifier

Poly magnifier with 4X magnification. 1 each L/R.

### MB015 Video Camera Arm

Easy positioning of camera over full length of grossing station via an aluminum slide extrusion. The overhead positioning is included with an articulating arm to support most camera systems. Video camera not included.

### MB016 Halogen Task Lighting - 2 each


Two halogen lighting system, 20 watts centered over the dissection table area. The overhead mount with 20-inch flexible arm allows for convenient and easy positioning.

### MB017 Camera Stand

Vertical camera stand with adjustable arm. Capable of supporting both digital and 35MM type cameras. Camera not included.

### MB019 Formalin Collection System

This system includes a collection funnel with perforated tissue screen piped to waste collection carboy located under the grossing workstation. Quick disconnects are included along with safety caps for safe transport to disposal or recycling. 2-1/2 gallon carboy with easy grip handling designed to

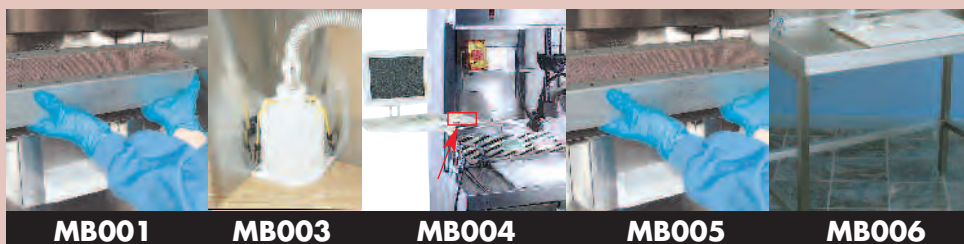

*Mopec offers an exceptionally wide variety of outstanding features. All of Mopec's standard workstations can be easily customized to meet the needs of your lab.*

## MB Series Optional Features

### Customized to Meet Any Laboratory Specifications

safely collect residual formalin for safe disposal or recycling. Formalin not included.

#### **MB020 Literature/Form Holder**

Lexan literature dispenser capable of holding 1" letter-size documents for easy handling and referral (three tier).

#### **MB026 Ventilated Trash Receptacle**

Removable ventilated trash container, 12 gallon capacity. Included with a clear flexible ventilation duct, which is tied into the grossing station's ventilation system.

#### **MB028 Formalin Dispensing System**

Dispensing system that is contained within the grossing station. Includes 5-1/2 gallon dispensing carboy, pump and backsplash mounted valve. Indicator lights for high/low formalin liquid levels in carboy. Formalin not included.

#### **MB029 Computer CPU Bracket**

Adjustable size mounting bracket for CPU, keyed to mounting plate with two pins for easy removal or installation. Computer not included.

#### **MB031 Removable Writing Ledge**

Add an additional writing ledge surface to your grossing station with the twist of a knob. Just slide it over the sink and you have an extended writing or working ledge that is easily removable.

#### **MB041 Combination Ledge and Pencil Drawer**

Convenient pull out writing ledge and utility drawer made of stainless steel construction.

#### **MB047 Pneumatic Foot Pedal**

Pneumatic foot operated water control. Foot operated toggle on/off switch. Hand operated mixing valve next to faucet.

#### **MB048 Stainless Steel Shelf**

Stainless steel shelf adjusts to accommodate storage needs.

#### **MB049 New LCD Flat Screen Monitor & Keyboard Adjustable Arm System**

Allows positioning of LCD display for greater visibility Up, down, forward and back. Full range of motion for maximum comfort level keyboard adjustable arm. Effortlessly position a keyboard in the most comfortable typing position, seated or standing use. Flat screen, monitor, keyboard and computer not included.

#### **MB050 Ventilated Trash Receptacle on Wheels**

Removable ventilated trash container, 12 gallon capacity with 3" swivel casters with brake mechanism. Included with a clear flexible ventilation duct, which is tied into the grossing station's ventilation system.

#### **BL800 Waste Disposal**

1 HP commercial grade disposal is air-cooled for continuous operation. Power switch activates solenoid to continuously feed water to the unit. Stainless steel internal parts prevent corrosion. Swivel impellers ensure the disposal will never jam.

*Mopec has developed a complete line of grossing workstations to fit your individual laboratory requirements. Each workstation can be customized with a variety of available options—all at an affordable price.*

