

CE900 Autopsy Table Pedestal



A hot and cold water faucet is provided with a single lever control and a vacuum breaker to safeguard the facility water supply. Shown with optional foot pedal.

The CE900 Autopsy Table is a pedestal style table of all stainless steel construction that incorporates all of the convenient features necessary to efficiently perform an autopsy. The CE900 is the shortest version of pedestal style autopsy tables and is ideal for pediatric institutions. This style also incorporates an integrally constructed sink into the table top as opposed to a sink wing. Unit comes standard equipped with Mopec's standard downdraft ventilation system.

A hydro-aspirator is conveniently located near the sink basin for ease in collection and/or disposal of bodily fluids. Electrical receptacles are provided on both sides of the table. A deck mounted sink spray with a high pressure hose allows complete rinsing of the body or table with the press of a single lever. The autopsy table has an integrally constructed large sink and is creased toward the center for complete drainage.

CE900 Standard Features

Recessed Table Top: The table top is fabricated of 14 gauge stainless steel with a #4 satin finish with built in large sink on one end.

Electrical: Two G.F.C.I. electrical receptacles are positioned on each side of the pedestal.

Hydro-Aspirator: Located on the end of the table near the sink basin. The aspirator has a reversing valve to clear the aspirator of any foreign matter. The aspirator is designed to drain into the sink for collection or disposal and is completely plumbed with control valve, vacuum breaker and tubing.

Faucet: Swing spout faucet has a vacuum breaker and single lever hot/cold handles.

Integral Construction: Integral sink to front panel

CE900 Autopsy Table Pedestal



CE900 shown with optional foot pedal


construction eliminates sharp edges where bacteria can accumulate. All welded seamless construction provides durability.

Sink Spray: Sink spray has chrome plated brass with a high pressure hose.

Pedestal: Table pedestal is fabricated of 16 gauge stainless steel with a #4 satin finish.

Grid Plate: Grid plates are recessed 1/2" and are supported on the table edge. The grid plates are perforated with 3/8" diameter holes on 2" centers. Four independent grid plates are perforated and removable.

Product No. CE900

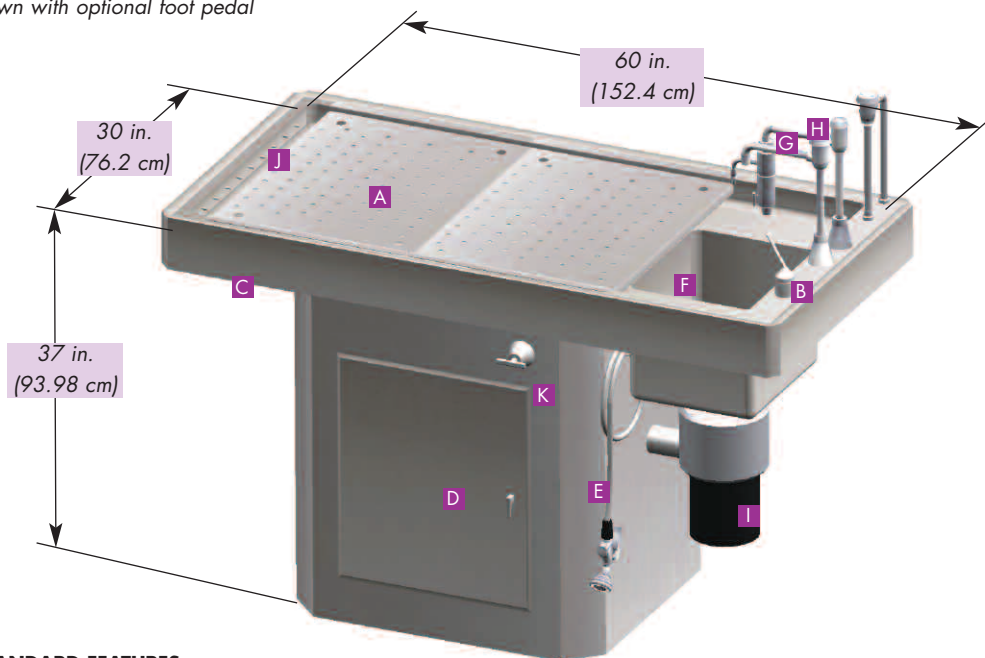

STANDARD FEATURES:

- A** Perforated Grid Plates (3 each)
- B** Hot/Cold Mixing Valve
- C** GFCI Outlet (2 each)
- D** Access Panel (2 each)
- E** Hand Spray Hose
- F** Sink Well
- G** Water Fixture with Swing Spout
- H** Hydro-Aspirator with Reverse Valve

OPTIONAL FEATURES:

- I** BL800 1 HP Disposal
- CO009 Perimeter Rinse Assy. (Not Shown)
- CO017 Scale Stand Table Mounted (Not Shown)
- J** CO023 Table Rinse Control Valve with Ten
- K** Nozzles (2 letters shown)