

OA190 Series Large Animal Cart

For Dissection/Transport



OA190 can easily transport large animals

This cart is a mobile unit that is excellent for transporting large animals throughout your facility. It is capable of carrying loads up to 1500 lbs, and is easily moved via a removable handle and a heavy-duty 6" diameter caster system featuring two rigid casters in back and two swivel casters in the front. A rear foot-operated brake brings the unit to a stop and keeps the cart in place while not in use.

Designed with a heavy-duty stainless steel top and large profile marine edge to help contain fluids, the table has a drain with a removable perforated basket strainer to collect debris.

Mopec offers design flexibility for the table top. Depending upon your needs, the cart can be rectangular or kidney shaped (OA110) allowing for easy access to large animals. A standard rectangular table top is 48" wide and 84" long.

Standard Features

Cart top: The cart top consists of dual 11 gauge stainless steel with a #4 satin finish. The top piece is fabricated with a 3-1/2 inch perimeter top edge raised at an angle of about 60°.

Cart perimeter: The center of the cart is creased toward a drain.

Construction: Utilizing TIG welding on all seams and mating sections, exposed welds are ground smooth and grained to match the workstation's #4 satin finish.


Handle: Removable 1-1/2" square stainless steel tubing.

Product No. OA190 — Kidney Shaped

Product No. OA191 — Rectangular

OA190 Series Large Animal Cart

For Dissection/Transport

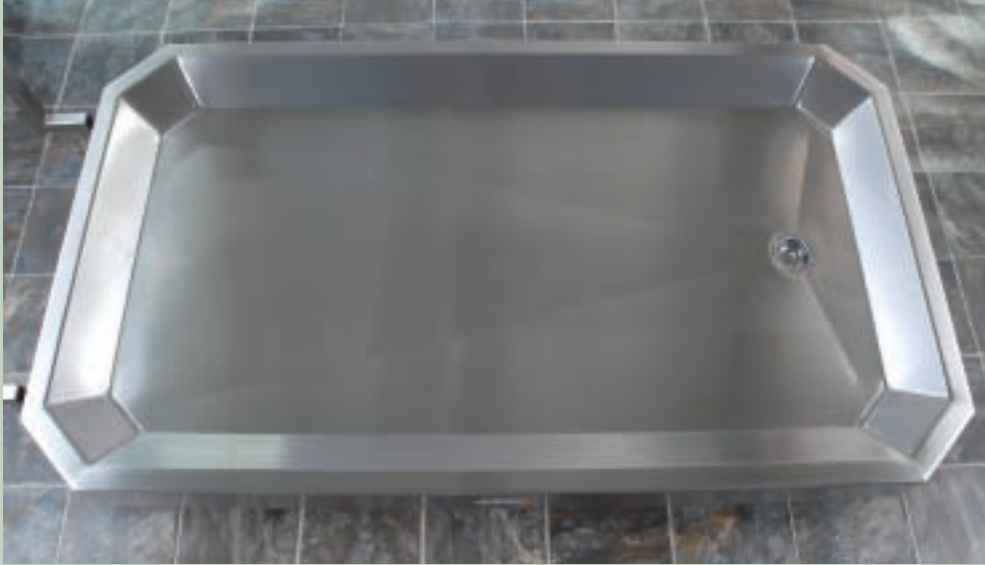

Top view of OA190

STANDARD FEATURES:

- A** Recessed table top (3-3/4" recess)
- B** Natural slope to the drain
- C** Drain hole with perforated basket strainer

