

OA500 Necropsy Table

Rotating Pedestal Mount Necropsy/Dissection Table



The OA500 Necropsy/Dissection Table is a pedestal design that incorporates all of the features necessary to efficiently complete a necropsy. Suited for smaller work areas, the OA500 necropsy table is shorter than conventional models and rotates allowing more versatile design applications. Mopec exclusively utilizes 16 stainless steel bearings for a smooth, accurate and virtually effortless 180 degree motion. Standard design features include Mopec's exclusive downdraft ventilation system that draws air into the center pedestal of the table. Also equipped with a pedestal anti-flood device which is an important feature since the building vent system is tied directly into the pedestal below the floor.



Table rotates



Table rotates

Check out
MOPEC'S STANDARD
FEATURES INCLUDE
Downdraft
Ventilation
Taking Care of Your Lungs

OA500 Necropsy Table

Rotating Pedestal Mount Necropsy/Dissection Table



Shown here without perforated grid plates

Standard Features

Table Top: The table is fabricated of 14 gauge type stainless steel with a #4 satin finish.

Electrical: An electrical receptacle is positioned on each side of the pedestal. (2 each)


Hydro-Aspirator: Located on the end of the table. The aspirator has a reversing valve to clear the aspirator of any foreign matter. The aspirator is designed to drain into the sink's recessed top for collection or disposal and is completely plumbed with control valve, vacuum breaker and tubing.

Faucet: Swing spout faucet has a vacuum breaker and a single lever hot/cold control mixing valve.

Pedestal: Table pedestal is fabricated of 16 gauge stainless steel with a #4 satin finish.

Perforated Grid Plate: Construction fabricated of 14 gauge stainless steel with a #4 satin finish.

Product No. OA500

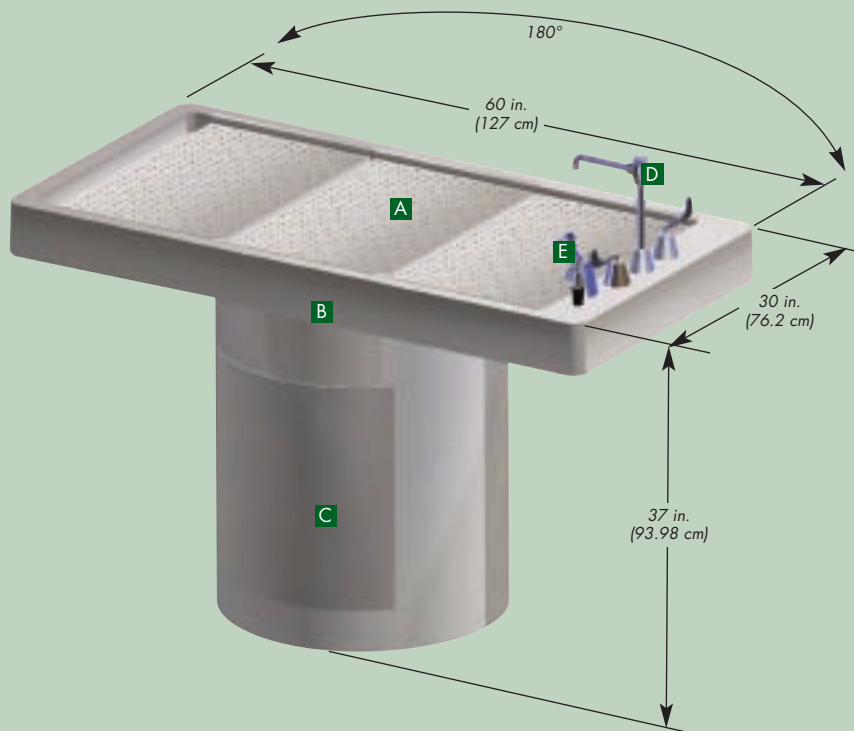

STANDARD FEATURES:

- A** Grid Plates
- B** Downdraft Cones
- C** Access Panel
- D** Water Fixture (H/C)

OPTIONAL FEATURES:

- E** OO010 Hand Spray Rinse
- OO025 Perimeter Rinse (Not Shown)
- OO026 Scale Stand Table Mounted (Not Shown)