


# OL100 Downdraft Ventilated Bedding Disposal (Dump) Station


Featuring Downdraft Ventilation



- Portable
- Safe
- Ergonomic

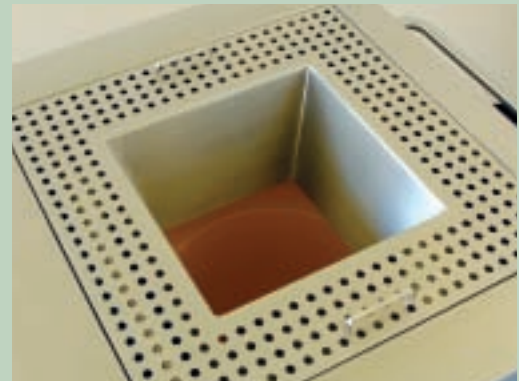

The OL100 has powerful downdraft ventilation

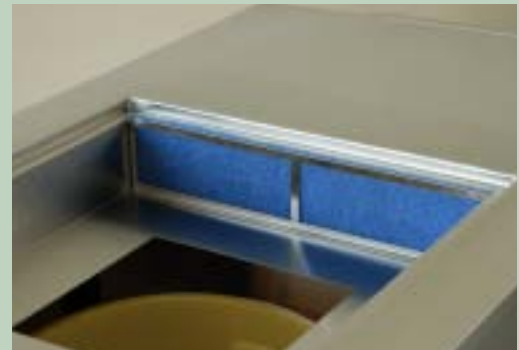

The OL100 features a two stage Hepa filtration system

Mopec's OL100 Downdraft Ventilated Bedding Disposal Station works with virtually any animal cage. Constructed of all stainless steel with an extra large (35" x 23") (88.9 cm x 58.42 cm) staging work surface. Access doors can be locked for safety. A self contained dust control fume evacuation system combined with a powerful Hepa filtration system offers features necessary for a clean, safe and efficient environment. Downdraft vent cone 15" x 15" (38.1cm x 38.1 cm).

## Standard Features

Table Top: The table top is fabricated of 16 gauge stainless steel.

Handles: Stainless steel push handles on both sides of station.

Casters: Three inch casters, lockable and directional.

Downdraft System: Mopec's downdraft system is a unique system of ventilating harmful and noxious odors during cage cleaning. Air evacuation is accomplished through down draft ventilation cone and the dual filtration system.

Bucket: 20 Gallon portable braced bucket and stand enclosed in air chamber.

Powerful Hepa Filtration System: Dust control fume evacuation.



The OL100 is constructed of all stainless steel for easy cleaning

# **OL100 Downdraft Ventilated Bedding Disposal (Dump) Station**

**Featuring Downdraft Ventilation**



*Portable Dumping Station features a 20 gallon mobile bucket*



*Extra large Hepa filtration with full service access*

